To validate your warranty please fill out the form inside of this manual or register online at [www.LaserComb.net/warranty](http://www.LaserComb.net/warranty)

[www.HairMax.com](http://www.HairMax.com)

Patented in the USA. Patent Pending in 104 Countries.  
HairMax LaserComb is a Registered Trademark of Lexington International, LLC.

## USER MANUAL

Proudly Patented and Manufactured in the USA by  
**Lexington International, LLC**

**Warranty information  
card inside**

**To validate your warranty please register  
online at [www.lasercomb.net/warranty](http://www.lasercomb.net/warranty)**

## TABLE OF CONTENTS

Glossary.....	3
Introduction.....	4
Contraindications.....	5
Warnings.....	5
Precautions.....	6
Electrical Safety.....	7
Indications.....	7
Risks.....	8
Benefits.....	8
Usage Recommendation.....	9
Summary of Clinical Data.....	9
Information on Condition.....	10
Norwood Hamilton Classification System.....	10
HairMax LaserComb and Parts.....	11
Instructions for Use.....	12
Safety Label.....	14
Power Supply Specifications.....	14
Maintenance Instructions.....	15
Warranty.....	16
How to Monitor Device Activity.....	17
Device Failure Mode.....	17
Expected Life of the HairMax LaserComb.....	17
Troubleshooting.....	18

## GLOSSARY

**Crown** - The top or highest part of the head.

### **Fitzpatrick Skin Type Classification Guide:**

A system that classifies your skin type based on your complexion and responses to sun exposure:

Type 1: Highly sensitive, always burns, never tans.

Type 2: Very sun sensitive, burns easily, tans minimally.

Type 3: Sun sensitive skin, sometimes burns, slowly tans to light brown.

Type 4: Minimally sun sensitive, burns minimally, always tans to moderate brown.

Type 5: Sun insensitive skin, rarely burns, tans well.

Type 6: Sun insensitive, never burns, deeply pigmented.

**Hair Density** - The amount of hairs on the head. Generally, it is measured by counting the number of hairs found in one square inch of scalp.

**LASER** - Light Amplification by Stimulated Emission of Radiation is an optical source that emits photons in a coherent beam.

**Low Level Laser** - Also known as soft laser, cold laser or laser biostimulation. The HairMax LaserComb uses this type of laser.

**Mean Terminal Hair Density** - The average number of thick, pigmented hairs in the specific area.

**Non-consecutive days** - A day in between the treatment days.

**Norwood Hamilton Scale** - System that classifies the different ways in which baldness progresses in men.

**Placebo Device** - A device that looks and sounds like the HairMax LaserComb. The placebo device does not have a laser. This device was used in the study to test the effectiveness of the HairMax LaserComb.


**Terminal hairs** - thick, pigmented hairs

**Vellus hair** - fine, non-pigmented hair

# HairMax LaserComb®

The HairMax LaserComb® is a laser-based device that helps promote hair growth. The HairMax LaserComb provides the convenience of treatment in your own home.

The HairMax LaserComb works directly on your scalp by parting your hair with two rows of teeth. The teeth are aligned with a row of laser beams. The teeth part your hair forming a clear path for the light to reach the scalp as you use the HairMax LaserComb.

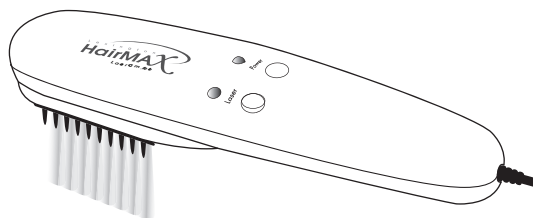This path is a critical feature to the success of the treatment. Without the teeth your hair would block the laser beams from reaching the scalp. Parting the hair provides direct access to the scalp.

Like brushing your teeth to maintain healthy teeth & gums, regular HairMax LaserComb usage is required. Three times a week on non-consecutive days for 10 to 15 minutes a session is all that is needed to achieve results over time. Over the next few months, you may notice a gradual increase in the number of hairs in a specific location (Hair Density). Most users notice continued improvement as they continue to use the HairMax LaserComb. If you stop treating your scalp with the HairMax LaserComb after you see hair growth, this improvement may gradually disappear.

## CONTRAINDICATION

If you are taking any photosensitizing medications or receiving any therapy that sensitizes the skin to light, DO NOT use the HairMax LaserComb. Please consult your physician if you have any other concerns.

## WARNINGS

- **WARNING:** Do not allow children to use or play with the HairMax LaserComb. Keep out of reach of children. It is not a toy.
- **WARNING:** Avoid getting the HairMax LaserComb wet. Contact with water may result in an electrical shock.
- **WARNING:** Never use the HairMax LaserComb if the cord or plug is damaged. Such use can cause the HairMax LaserComb to malfunction and produce an electrical shock. If the HairMax LaserComb is not working properly or if it has been dropped, damaged or gotten wet, contact Lexington International, LLC at 561-417-0200.
- **WARNING:** Avoid using the HairMax LaserComb near heated surfaces. Contacting a heated surface with the cord or plug can cause the product to malfunction and produce an electrical shock.
- **WARNING:** Only use the HairMax LaserComb according to the instructions provided. Any other use is not recommended.
- **WARNING:** Laser Light - Avoid Direct Eye Exposure. Do Not point laser or allow laser light to be directed or reflected toward other people. Do Not stare directly at the laser light or its reflection off a mirror-like surface. Doing so can cause temporary irritation to your eyes. Prolonged direct exposure can harm your eyes.

## PRECAUTIONS

- **CAUTION:** Never leave the HairMax LaserComb unattended when plugged into a power outlet. An electrical malfunction can result in a fire.
- **CAUTION:** Do Not put strain on the cord when untangling or straightening it. Doing so can break the cord causing it not to work.
- **CAUTION:** Never use magnifying optical instruments with the HairMax LaserComb. Their use with this product can cause temporary irritation to the eyes.
- **CAUTION:** Never use acetone or other solvents on any part of the HairMax LaserComb. Use of solvents on the window can block the laser light from coming through.
- **CAUTION:** Not recommended for use by men with Norwood Hamilton Classification I, II, Va, VI, or VII (See drawings on page 10 of the User Manual) because they are not likely to experience hair growth from using the HairMax LaserComb.
- **CAUTION:** There is currently no validated scientific evidence supporting the safety and effectiveness of the HairMax LaserComb on women.
- **CAUTION:** The safety or effectiveness of the HairMax LaserComb in men with Fitzpatrick Skin Types V and VI (see the description on these skin types on page 3 of the user manual) has not been evaluated.

**NOTICE:** If any malfunction occurs, do not attempt to repair it yourself. This appliance has no user-serviceable parts and is not to be opened or modified by the user. Use the HairMax LaserComb only for its intended use as described in this manual. Contact Lexington International LLC for repairs at 561.417.0200.

## ELECTRICAL SAFETY

Basic safety precautions should be observed when using this, or any other electrical device, especially when children are present.

**⚠ DANGER:** Keep away from water. Contact with water may result in an electrical shock.

To reduce the risk of hazard by electric shock:

Always unplug the HairMax LaserComb from the power outlet immediately after using.

Do not use while bathing.

Do not place or store the HairMax LaserComb where it can fall or be pulled into a tub or sink.

Do not expose the HairMax LaserComb to running water, or drop into water or any other liquid. If the HairMax LaserComb falls into water, unplug it immediately from the power outlet.

Do NOT reach into or touch the water.

If using the HairMax LaserComb near water, the use of an electrical outlet protected by a Ground Fault Circuit Interrupter is recommended.

## INDICATIONS

The HairMax LaserComb is indicated to promote hair growth in males with androgenetic alopecia who have Norwood Hamilton Classifications of IIa to V and Fitzpatrick Skin Types I to IV.

## RISKS

Harm to Eyes – Do Not shine the light from the HairMax LaserComb directly into your eyes. **Prolonged** direct exposure to your eyes can cause harm to your eyes.

Scalp Itching – Less than 5% of the subjects in the clinical trial of the HairMax LaserComb reported mild scalp itch during the usage of the HairMax LaserComb. If you experience more than mild itching consult a physician.

Tingling – Less than 4% of the subjects in the clinical trial reported a mild and temporary tingling feeling. This lasted for a short period of time after applying the treatment. If the problem persists consult a physician.

Increased blood flow to the scalp – Less than 1% of the subjects in the clinical trial reported a mild increase in blood flow (**also known as vasodilatation**) to the scalp. This lasted for a short period of time after applying the treatment. If the problem persists consult a physician.

## BENEFITS

The HairMax LaserComb is a home-use Low-Level Laser treatment to promote hair growth in men. The HairMax LaserComb uses laser light to treat and stimulate your scalp.

Individual results may vary. Your hair grows at a slow pace. Increases in the number of hairs in one specific location (**hair density**) are difficult to notice and are not immediately visible. A typical user will notice some increase in hair density within the first few weeks. Over time you will begin to re-grow hair.

Proper use according to the instructions leads to positive results. If you stop treating your scalp with the HairMax LaserComb after you see hair growth, this improvement may gradually disappear.

## USAGE RECOMMENDATION

Use the HairMax LaserComb 3 times per week with at least a day between each treatment (**non-consecutive days**). (For example Monday, Wednesday, and Friday.) The total time of each treatment is 10-15 minutes. More frequent or longer usage will not increase results. Maximum results will not be seen if you do not follow the recommended treatment schedule.

Like brushing your teeth to maintain healthy teeth and gums ongoing HairMax LaserComb usage is needed to maintain healthy hair growth. If you stop using the HairMax LaserComb, any benefit may gradually disappear.

### Summary of Clinical Data


A study was conducted at 4 different locations around the country. All participants in the study (ages 30-60) were males, had hereditary hair loss, which is called **Androgenetic Alopecia**, and light to medium skin on the head. A total of 123 men participated in the study. Subjects were divided into two groups. One group used the HairMax LaserComb. The second group used a **placebo device** that looked and sounded like the HairMax LaserComb but did not include a laser. The first group included 82 participants. The second group included 41 participants.

Subjects performed treatments 3 times per week on non-consecutive days for 6 months. Hair growth was monitored during the study. The number of thick, pigmented (terminal) hairs in a specific area were compared from the beginning to the end of the study. 67 of 72 (93%) men using the HairMax LaserComb had some increase in the number of terminal hairs in that area. The average number of terminal hairs per square centimeter (**mean terminal hair density**) measured in this study increased by 19 over a 6 month period. The area measured was roughly the size of a dime. No subjects experienced serious adverse events. The number and types of adverse events were similar in both groups.

## Information on Condition


Androgenetic Alopecia (AGA) is a condition involving both genetic (inherited) and hormonal factors. AGA in men usually begins at the hairline. The hairline gradually recedes to form an "M" shape. The hair becomes finer and shorter. The hair at the top of your head (also called the **crown**) begins to thin. Over time the top of the hairline meets the thinned crown, leaving a pattern of hair around the sides of your head. A system called the Norwood Hamilton classification system is used to determine the extent of hair loss. The following chart will assist in determining your stage of hair loss.

### Norwood Hamilton classification system

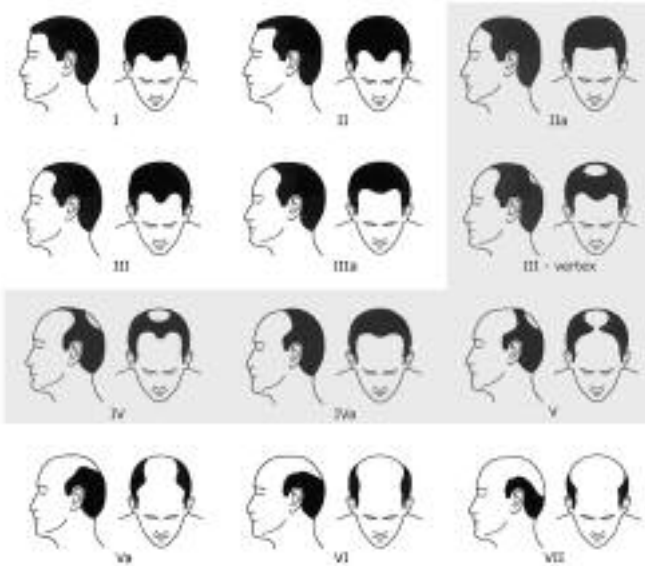An estimated 50 million men in the United States experience AGA. A person's tendency to lose hair can be inherited from either side of a person's family or from both parents. Men of all races and ethnicities lose their hair. By age 40, nearly forty percent of men have visible signs of AGA. By age 50, fifty percent of men show signs of the condition.

# HAIRMAX LASERCOMB AND PARTS

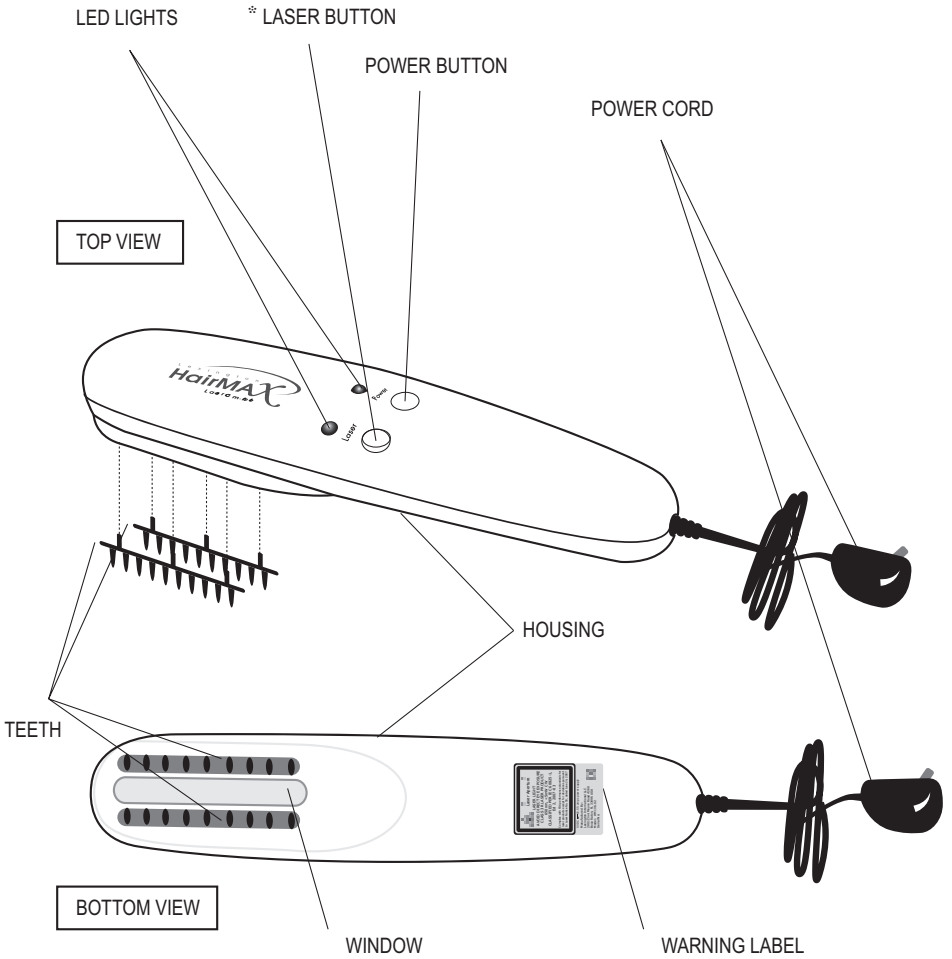\* Laser Button - Caution - use of controls or adjustments or performance of procedures other than those specified herein may result in hazardous radiation exposure.

## HAIRMAX LASERCOMB INSTRUCTIONS FOR USE




The HairMax LaserComb produces results if used correctly. It is manufactured using precision components and should last 10-15 years (see page 17 for details). Please take time to read this entire instruction manual before operating the HairMax LaserComb.

Use the HairMax LaserComb for 10 to 15 minutes per treatment, three times a week on non-consecutive days.

1. Wash your hair. Hair can be wet, towel dried or dry. The HairMax LaserComb should be used before applying any other topical treatments, gels, sprays or concealing products.
2. Insert the power attachment into the handle end of the HairMax LaserComb (See Illustration 2).
3. Plug the power adapter into any convenient wall electrical outlet (See Illustration 3).
4. Choose whether you want to use the beep tone.

**Caution** - use of controls or adjustments or performance of procedures other than those specified herein may result in hazardous radiation exposure.

To turn on the HairMax LaserComb and the beeping tone, push the 'Power' button, followed by the 'Laser' button (See Illustration 4). An audible beeping tone will be heard every 4 seconds. This tone helps you determine the time period between movements.

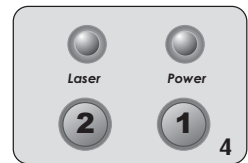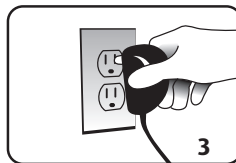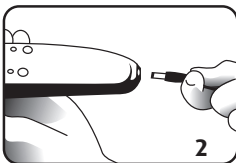
If you prefer to use the HairMax LaserComb without the beep tone, press the 'Laser' button a second time. You will hear a beep response and the audible tone will stop. To re-activate the tone, simply press the 'Laser' button again.

5. Place the tips of the HairMax LaserComb teeth flat on your scalp (See Illustration 5).
6. Move the HairMax LaserComb slowly by placing the HairMax LaserComb on a spot, leaving it there for 4 seconds then gliding it 1/2 inch to the next spot. Do not lift the HairMax

LaserComb off your head. Move the HairMax LaserComb against the direction of the hair growth - front to back on the first pass covering your entire head (see illustration 6).

**IMPORTANT:** Do not lift the HairMax LaserComb off of your head when moving every 4 seconds. The HairMax LaserComb features a pair of comb teeth that part your hair so that the light can reach the treatment areas. Lifting the HairMax LaserComb prevents your hair from being parted.

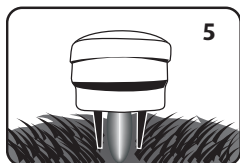7. Move the HairMax LaserComb slowly from one side of your head to the other side of your head on the second pass, covering your entire head (See Illustration 7).

**TIP:** On long hair it is best to use HairMax LaserComb as follows:

Hold the HairMax LaserComb in one hand and a normal comb in the other.

Place the comb upright and 1/2 inch in front of the HairMax LaserComb in the direction of the movement.

Lift up the hair with the comb as you move both the comb and the HairMax LaserComb progressively over your scalp (See Illustration 8)

The use of the comb helps move most of the hair out of the way of the HairMax LaserComb so that the laser beams can reach your scalp.

8. After treatment, turn off the power button, remove the wall plug.
9. Clean the window and comb teeth according to Maintenance Instructions found on page 15.
10. Return the HairMax LaserComb and power supply to the case for storage.

## Please Note:

The total time required for each treatment with the HairMax LaserComb is 10 to 15 minutes. It is important to move the HairMax LaserComb very slowly. You are bathing your scalp with laser light. If you move it too fast, you will not see the full benefit of the treatment.

### The following laser safety label can be found on this product:

**Wavelength:** 655 nm  
visible red light

This Laser Product is designated as Class 3R during all procedures of operation.

**Note:** HairMax LaserComb has no user-serviceable parts.

**Location:** On the underside of the handle.

Accessible  
Laser  
parameters

Wavelength	655 nm
Laser Power for classification	<5mW CW
Beam Diameter	<5 mm
Divergence	57 mrad
Number of lasers	9

Please visit <http://www.hairmax.com/QualitySafety.aspx>  
for more information or call 561.417.0200

## POWER SUPPLY SPECIFICATIONS

Output from transformer: 5V, minimum 500mA. Plug: 2.1mm x 5.5mm.

The power supply for the HairMax LaserComb can be used in the United States and internationally. This power supply will adapt to all power supplies from 110 volt/60 cycle (US) to 240 volt/50 cycle. If outlets in your country have different plugs than the power supply that comes with the HairMax LaserComb, use an inexpensive plug adapter not supplied with the device.

## MAINTENANCE INSTRUCTIONS

Always unplug the HairMax LaserComb from the power outlet while cleaning it.


DO NOT use acetone or any other solvents on any part of the HairMax LaserComb. Use of acetone or other solvents to clean the window can immediately block most of the laser light from coming through.

To clean the laser window (from where the laser light shines through), gently clean with a soft, slightly moistened cloth.

### Cleaning comb teeth

To keep the teeth of the combs clean, gently remove any hairs or dust with a slightly moistened towel.

For a more thorough cleaning, gently pull out comb teeth. Teeth can be washed using soap and water. Dry combs thoroughly, and replace them using a little gentle manipulation.

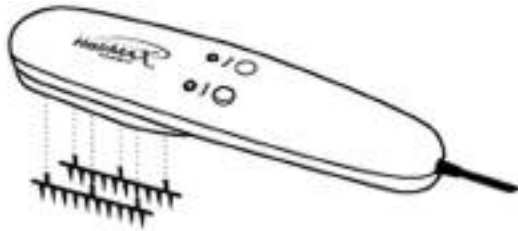Comb teeth should be kept clean but do not need to be replaced. To order replacement combs, please contact our office at 561.417.0200.

If the cord becomes twisted, unplug and straighten it prior to use.

Do not put strain on cord if you need to untangle or straighten it.

Only use directions, cleaning methods and operating procedures specified in this manual.

## **WARRANTY**

The HairMax LaserComb is warranted by Lexington International LLC ("Manufacturer") against defects in materials or workmanship for 24 months from date of purchase. If this product is determined to be defective, Lexington will repair or replace it, at its option, at no charge during the warranty period.

After the warranty period, the user must pay for all parts and labor charges. If the product is repaired or replaced, either within or outside the warranty period, Lexington reserves the right to provide equivalent parts in different colors to your original product.

This warranty does not cover damages caused by acts of God, misuse, negligence, accident, modification of, or to any part of the product. This warranty does not cover damage due to improper operation or maintenance, connection to improper voltage supply, or attempted repair by anyone other than Lexington. This warranty is void if product is tampered with.

To obtain warranty service, you must deliver the product freight prepaid, in either its original packaging or packaging affording an equal degree of protection, together with your name and address, to:

**Lexington International LLC  
777 Yamato Road, Suite 105  
Boca Raton, FL 33431, USA**

Limitation of Liability; Indemnity: Except as expressly set forth above, Lexington International LLC shall not be liable for any special, indirect, or consequential damages that are directly or indirectly related to the use of, or the inability to use, the HairMax LaserComb, even if Lexington International LLC has been advised of the possibility of such damages whether in an action of contract, negligence or other tortious action. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation or exclusion may not apply to you.

Governing Law: The usage of the HairMax LaserComb and the interpretation of its terms shall be governed by and construed in accordance with the laws of the State of Florida applicable to agreements negotiated, executed and performed wholly and solely within the State of Florida. The parties agree that the State of Florida's conflicts of laws rules shall not apply to this Agreement or any transaction performed by or under this Agreement.

Jurisdiction and Venue: The parties irrevocably submit and consent to the sole and exclusive jurisdiction and venue of the Florida state courts in and for Palm Beach County, Florida, U.S.A. and the Federal Courts in and for the Southern District of Florida, U.S.A. The parties agree not to raise the defense of forum non conveniens.

Authorized Representative in Europe (Regulatory Affairs Only): Emergo Europe, Molenstraat 15, 2513 BH The Hague, The Netherlands. Tel: (31) (0) 70 345-8570 Fax: (31) (0) 70-345-7299

**Congratulations on becoming a HairMax LaserComb user!**  
*To validate your warranty please register online at*  
*[www.lasercomb.net/warranty](http://www.lasercomb.net/warranty) or call 561.417.0200*

## **How to Monitor Device Activity**

Your HairMax LaserComb is inspected and tested during production. When working properly, both the green and the red lights above the buttons turn on. A red laser light shines through the window at the bottom of the device. Do NOT look directly at the laser light. Shine the laser light at the palm of your hand to see red laser light “speckles” appear. You will also hear a beep tone.

Do not operate the laser if the unit is defective or damaged. If one or both of the indicator lights is not lit, the laser light does not appear, or there is no beep tone, contact Lexington International LLC for repair or replacement.

## **Device Failure Mode**

The HairMax LaserComb is designed to last 10-15 years of continuous use. The laser light does not fade or diminish over time. At the end of the laser’s life or if it fails for any other reason no red light will shine out of the device. Do not look directly at the laser light. Direct the laser light at the palm of your hand to see red laser light “speckles” appear. If no light appears, contact Lexington International LLC for repair or replacement.

## **Expected Life of the HairMax LaserComb**

The HairMax LaserComb is a medical device manufactured under strict quality controls and validation through GMP (Good Manufacturing Practices). The HairMax LaserComb contains a high quality laser diode designed to last at least 15,000 hours. Based on that rating and the durability of the components being used, the expected life of the HairMax LaserComb is 10-15 years of regular use.

## TROUBLESHOOTING

### *My HairMax LaserComb will not turn on*

- Insert plug into the HairMax LaserComb. Press the power button then the laser button.
- Try another wall outlet

If the HairMax LaserComb still does not work contact Lexington International LLC for repair or replacement.

### *My HairMax LaserComb is turning off during use*

- Try another wall outlet

Contact Lexington International LLC for a replacement power adapter

### *I hear a loud rattle inside my HairMax LaserComb when I move it*

- Hold down the buttons while shaking the HairMax LaserComb

If you still hear a rattle contact Lexington International LLC for repair or replacement.

### *My HairMax LaserComb beeps irregularly or does not beep*

- Unplug HairMax LaserComb. Plug in adapter and press the power button then the laser button.

If problem persists, contact Lexington International LLC for repair or replacement.

### *My HairMax LaserComb beeps but I do not see the laser*

Contact Lexington International LLC for repair or replacement.

**Contact Information:** Lexington International, LLC  
777 Yamato Road, Suite 105  
Boca Raton, FL 33431, Tel. 561.417.0200  
service@lasercomb.net, www.lasercomb.net

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# Warranty Card

The HairMax LaserComb® has a 2 year warranty against defects in materials or workmanship from Lexington International, LLC. Items that are damaged or abused are not covered under the warranty and will be charged at their regular repair rates. Lexington may extend this warranty period under certain product failures. Please refer to the User Manual for complete warranty terms. We guarantee your satisfaction.

**To validate your warranty please return this completed warranty card or register online at: [www.hairmax.com/warranty](http://www.hairmax.com/warranty)**

**Name:** \_\_\_\_\_

**Serial Number:** \_\_\_\_\_

**Date of Purchase:** \_\_\_\_\_

**Address:** \_\_\_\_\_

\_\_\_\_\_

**City:** \_\_\_\_\_ **State:** \_\_\_\_\_

**ZIP/Postal Code:** \_\_\_\_\_ **Country:** \_\_\_\_\_

**Email:** \_\_\_\_\_

**Where did you purchase this product (Reseller Name)?**

\_\_\_\_\_

**What factors influenced your purchase of this product?**

\_\_\_\_\_

\_\_\_\_\_

\_\_\_\_\_

*\*Please fill out this form, fold it, staple or tape it at the designated area, place postage stamp and mail it back to us.*

or staple  
tape

fold here

---

---

---

---

**Place  
Stamp  
Here**

**Lexington International, LLC  
2901 Clint Moore Rd., #260  
Boca Raton, FL 33496  
USA**

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Copyright 2007 HairMax LaserComb, LLC.